

# BASIC INFORMATION

Straumann® Modular Cassette





# CONTENT

<b>1. STRAUMANN® MODULAR CASSETTE – SYSTEM OVERVIEW</b>	<b>2</b>
1.1 Straumann® Modular Cassette, A Module	4
1.2 Straumann® Modular Cassette, B Module	8
1.3 Straumann® Modular Cassette, C Module	15
1.4 Lids	17
<b>2. STEP-BY-STEP REPROCESSING</b>	<b>18</b>
2.1 Point of use	18
2.2 Containment and transportation	18
2.3 Pretreatment	19
2.4 Cleaning	21
2.5 Inspection	24
2.6 Packaging	25
2.7 Sterilization	26
<b>3. FURTHER INFORMATION</b>	<b>27</b>
<b>4. PRODUCT REFERENCE LIST</b>	<b>28</b>

# 1. STRAUMANN® MODULAR CASSETTE – SYSTEM OVERVIEW

The Straumann® Modular Cassette is used for the secure storage and reprocessing of surgical and auxiliary instruments of the Straumann® Dental Implant System. The Straumann® Modular Cassette works with any Straumann® implant line (e.g. SP, BLT™, BLX™), including with the Straumann® Guided Surgery workflow. The system consists of three modules named A, B and C.

- The A Module stores tools that can be shared among different implant lines. Removable trays provide dedicated spaces to store instruments.
- The B Module stores tools for a specific implant line. Removable trays are dedicated to an implant line workflow.
- The C Modules have embedded inserts to store additional instruments without removable trays.

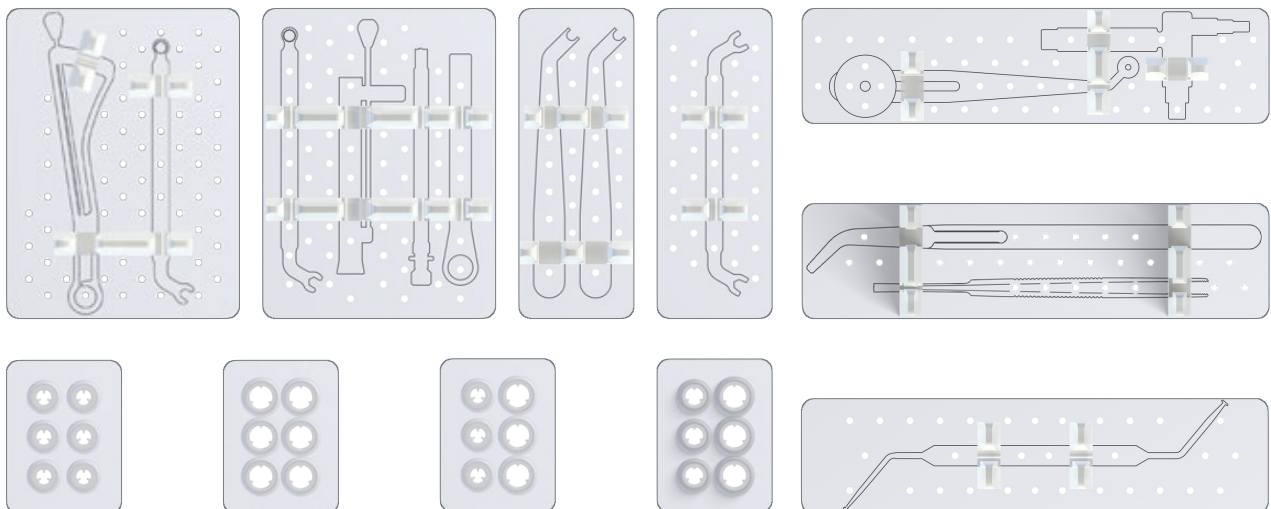
As the modules are not intended to maintain sterility on their own, they should be enclosed in a steam sterilization packaging system and sterilized.

Each module comes with a writable blank label plate to help differentiate each set. These label plates can be pushed out from the hole in the back for cleaning or replacement. Laser-marked label plates are also offered as spare label plates (Art. No. 041.775).


## Straumann® Modular Cassette, A Module (041.761)




## A Module trays




**Straumann® Modular Cassette, B Module (041.776)**

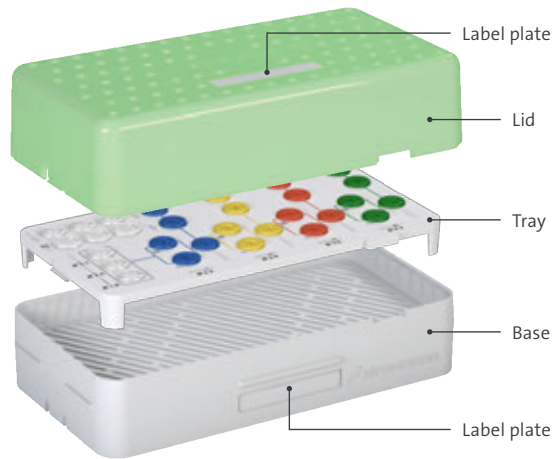


**B Module trays**

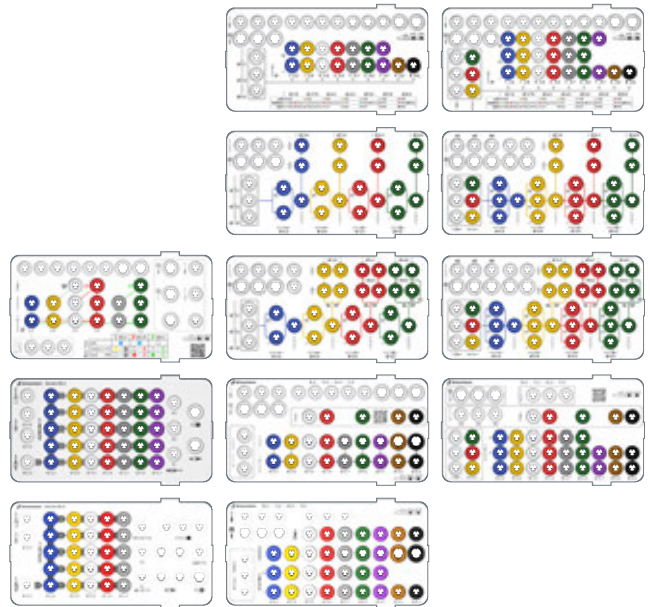


**Straumann® Modular Cassette C Module, Auxiliary Box (041.783)**

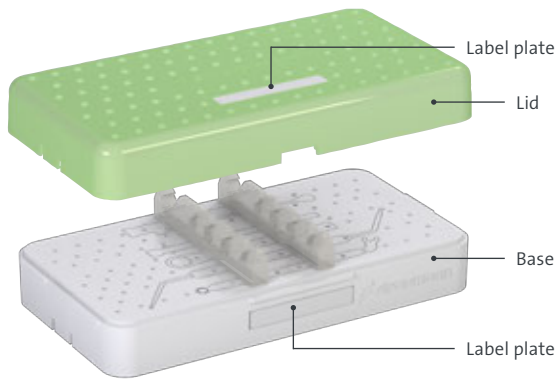

**Straumann® Modular Cassette C Module, Guided surgery (041.772)**

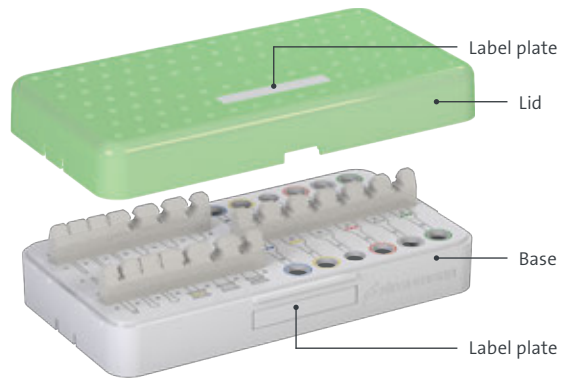


**Straumann® Modular Cassette C Module, Spare Tools (041.773)**

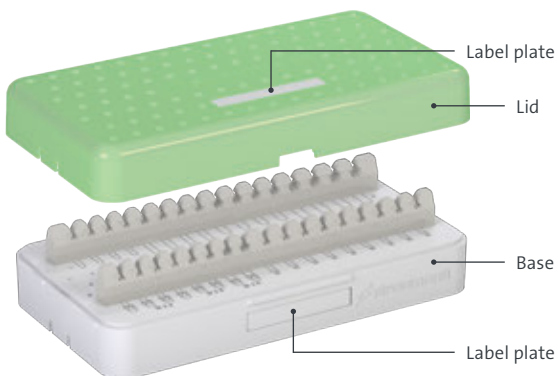


## 1.1 STRAUMANN® MODULAR CASSETTE, A MODULE

The A Module (Art. No. 041.761) mainly stores surgical tools that can be shared among different implant lines. Users can set up the A Module according to their needs by changing the removable trays inside the A Module. Refer to the *Straumann® Modular Cassette Selection Guide (702824/en)* to select the correct cassette set-up for your clinic.

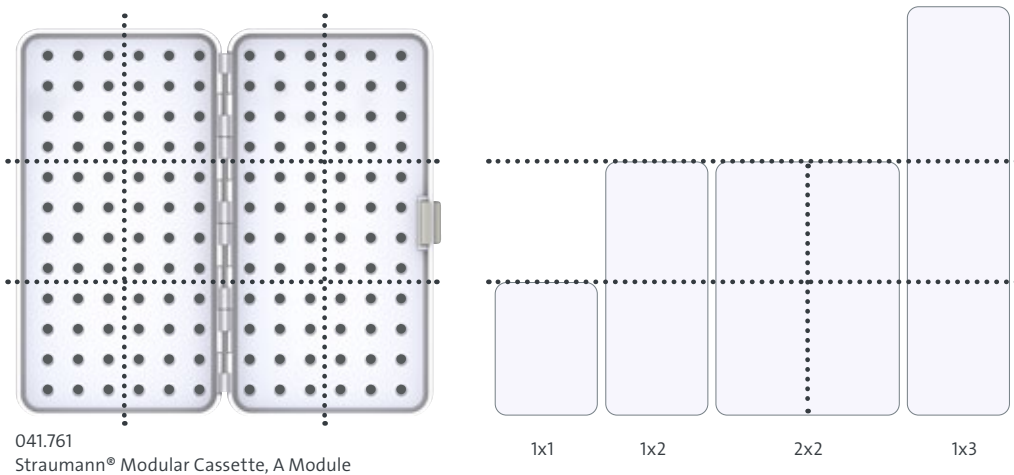The A Module box is a container for the trays. The grid allows the box to be set up as needed. There are four different tray sizes. Maximum loading capacity per box side is 2 by 3 of the 1x1 trays.



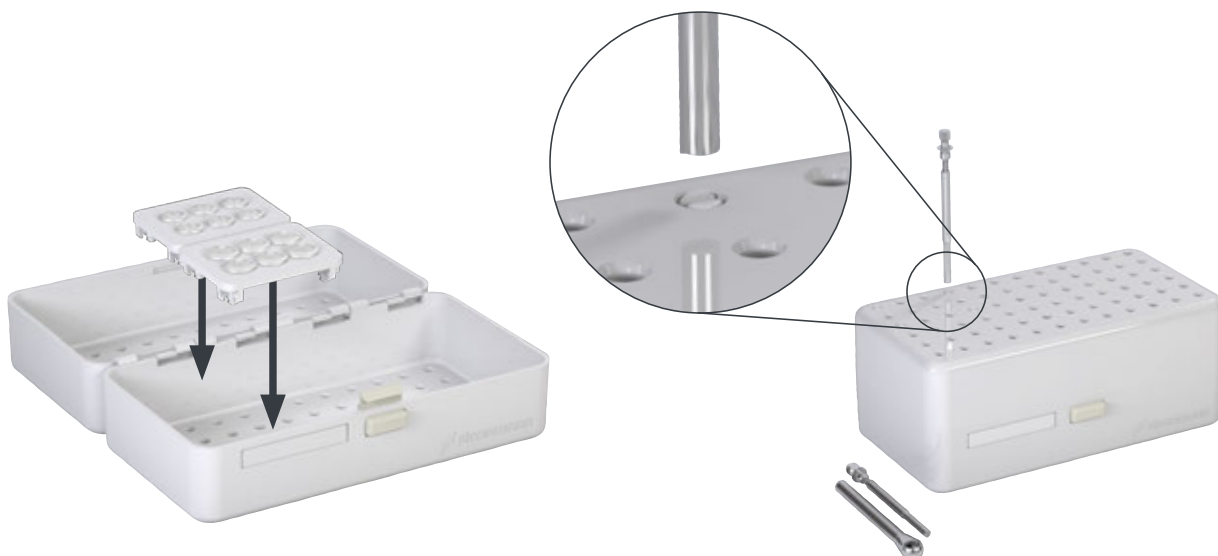To insert a tray, simply click the tray into the interior holes of the A Module. Physical feedback and an audible click will confirm secure placement of the tray. To remove, use a blunt instrument (e.g. the ratchet bolt from the disassembled ratchet, Art No. 046.119) to push the trays out from the back.






There are six different trays for long surgical instruments, each with silhouettes for easy identification of the different instruments held.



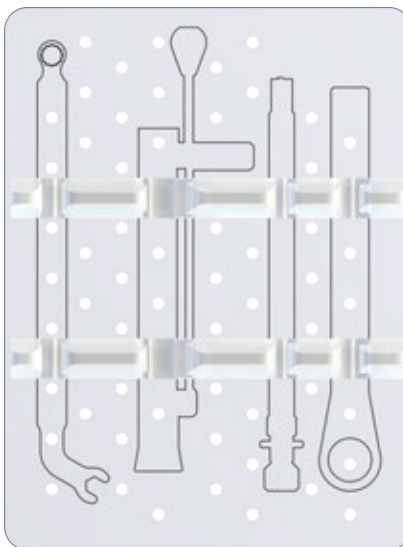
041.791  
Ratchet with Torque Control Tray




041.766  
Ratchet Tray



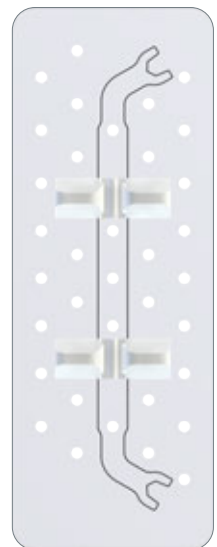
041.767  
Release Aid Tray



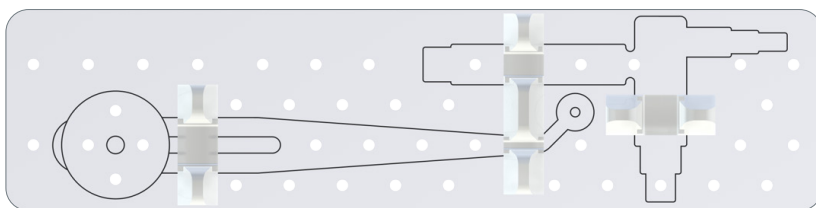
041.768  
Stop Key Tray



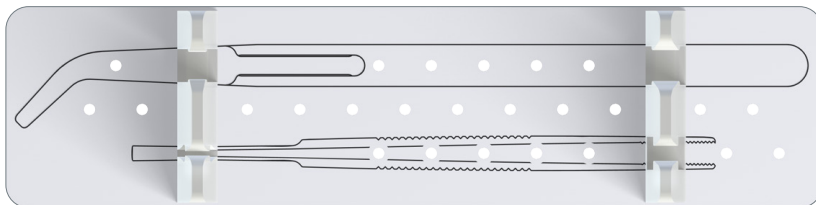041.769  
Planning Aid Tray



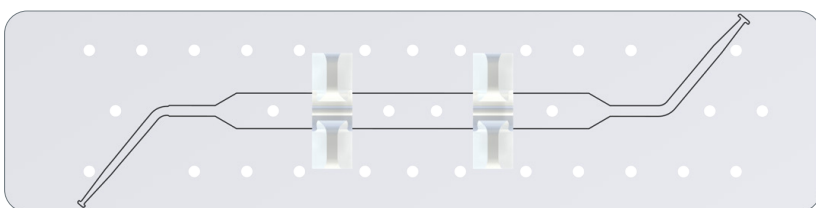041.770  
Tweezer Tray



041.771  
Implant Depth Gauge Tray



There are four different grommet trays. The grommet tray can store shorter tools such as

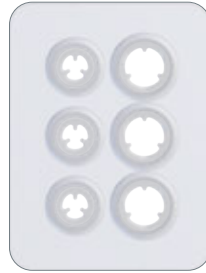
- SCS screwdriver
- Angulated screwdriver
- Implant adapter
- Bone profiler
- 48-hour explantation device.




041.762  
6 small



041.763  
6 large



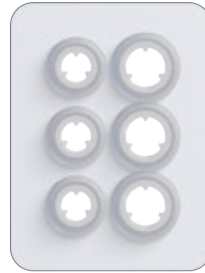041.764  
3 small, 3 large




041.765  
3 medium, 3 large

Users are free to assemble the A Module with any desired arrangement of tray combinations inside the module. Do not connect trays and label plates outside the A Module other than on the dedicated space for the label plate.

Verify that instruments do not touch each other when closing the A Module.

## Example of A Module setup

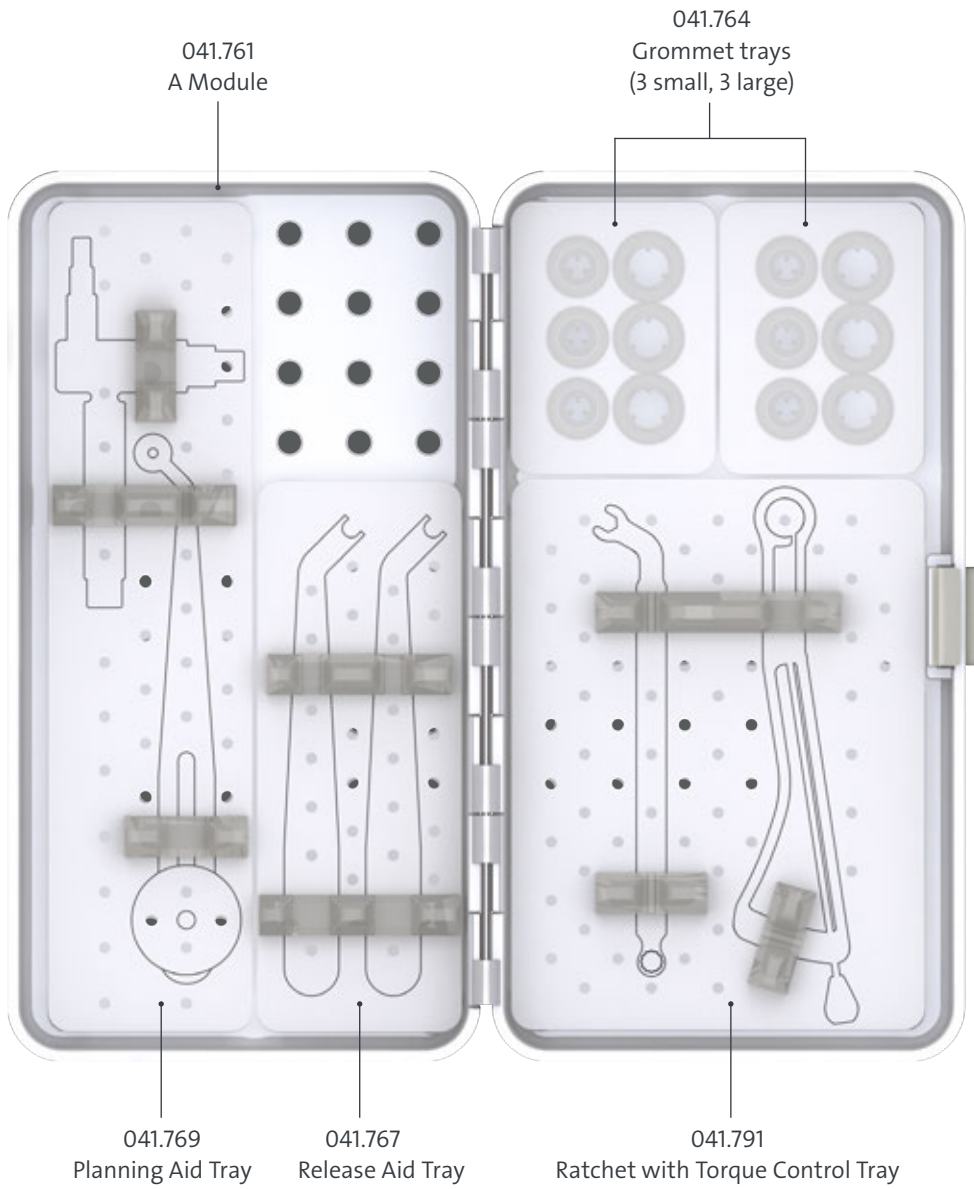


## 1.2 STRAUMANN® MODULAR CASSETTE, B MODULE


The B Module features different workflow trays that store cutting tools for specific implant lines. The B Module should be used together with an A Module to complete the tools needed for an implant surgery.

External slots on the side of the module help align each module on the surgical table. Refer to the *Straumann® Modular Cassette Selection Guide (702824/en)* to select the correct cassette set-up for your clinic.


To click the B Module trays in place, align the internal slots and gently push the tray into the base by hand. Push down on the connection point for trays. Physical feedback and an audible click will confirm secure placement of the tray. To remove the tray, gently pull the tray out by hand, holding the alignment tabs.
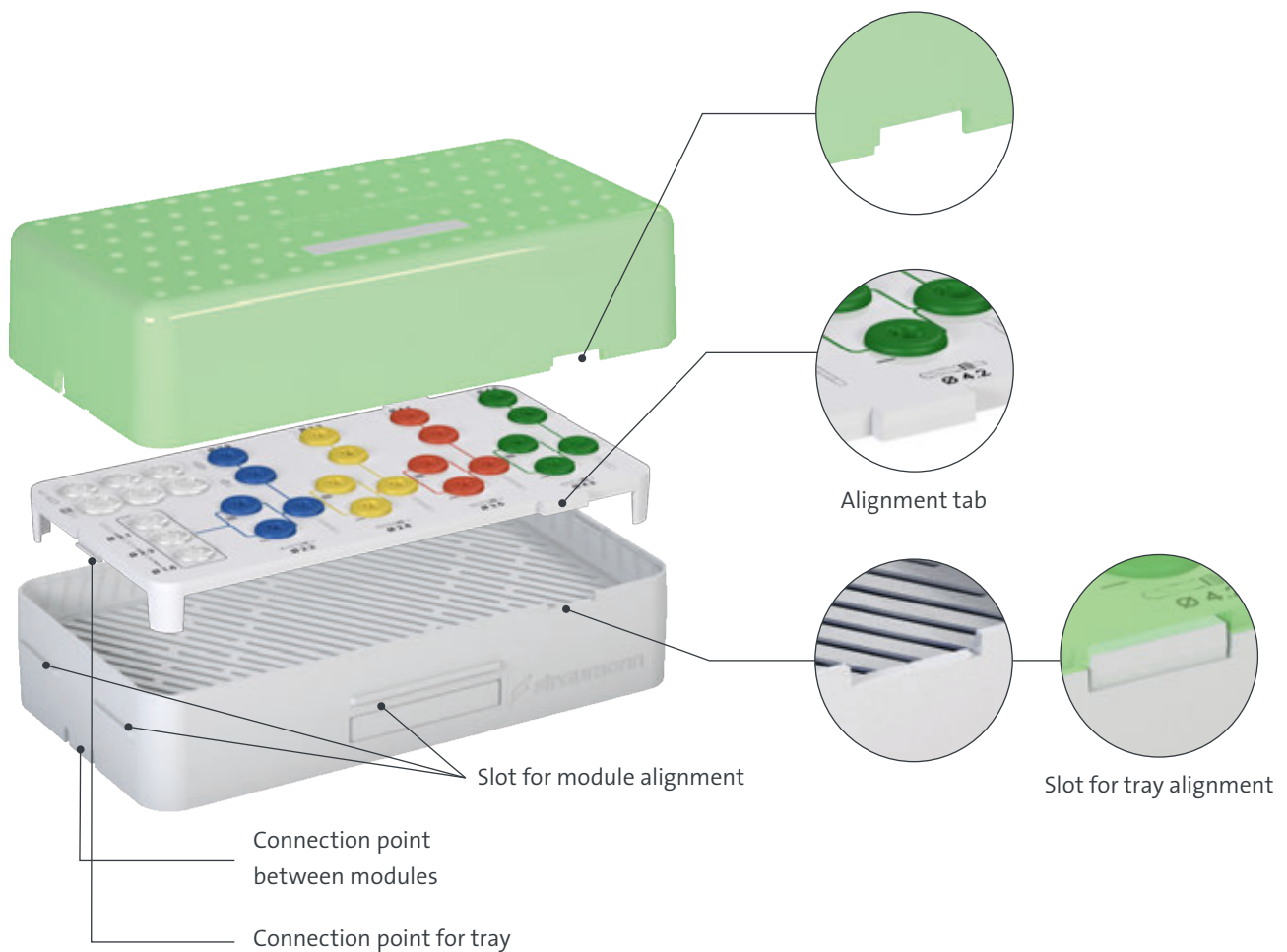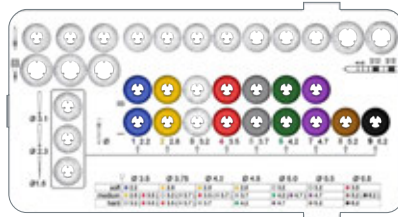
### Straumann® Modular Cassette, B Module (041.776)






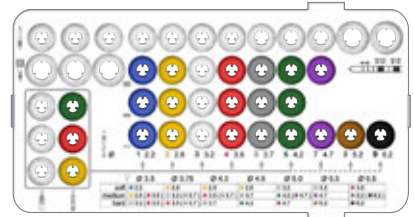
041.777 Fully Tapered Tray



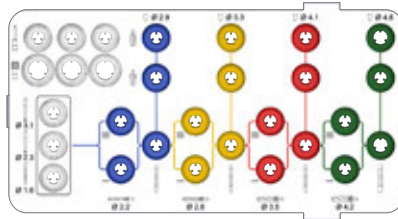
041.780 Guided Fully Tapered Tray



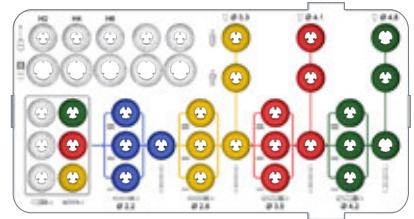
041.778 Basic Tray



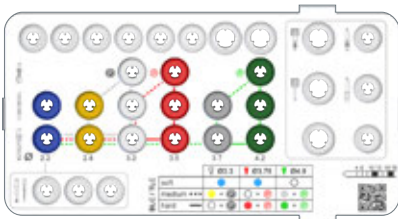
041.781 Guided Basic Tray



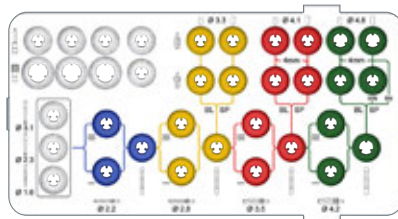
041.785 BLC, TLC Tray



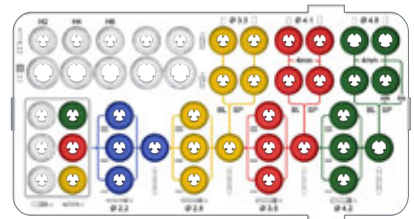041.779 Parallel Walled Tray



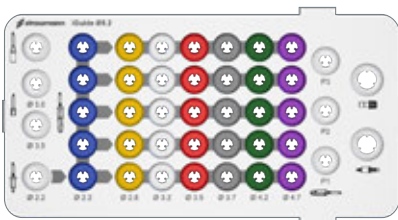041.782 Guided Parallel Walled Tray



041.789 Straumann iGuide™ B-Tray, Ø 5.2



041.787 BLC, TLC, BLX, TLX Tray



041.788 Guided BLC, TLC, BLX, TLX Tray



041.811 Straumann iGuide™ B-Tray, Ø 3.9




041.814 BLC, TLC, BLX, TLX 2 lines of drills Short & Long Tray

**1.2.1 Straumann iEXCEL™ Tray (BLX, TLX, BLC, TLC), Guided (041.787 & 041.814 & 041.788) and Straumann iGuide™ tray (041.789 & 041.811)**

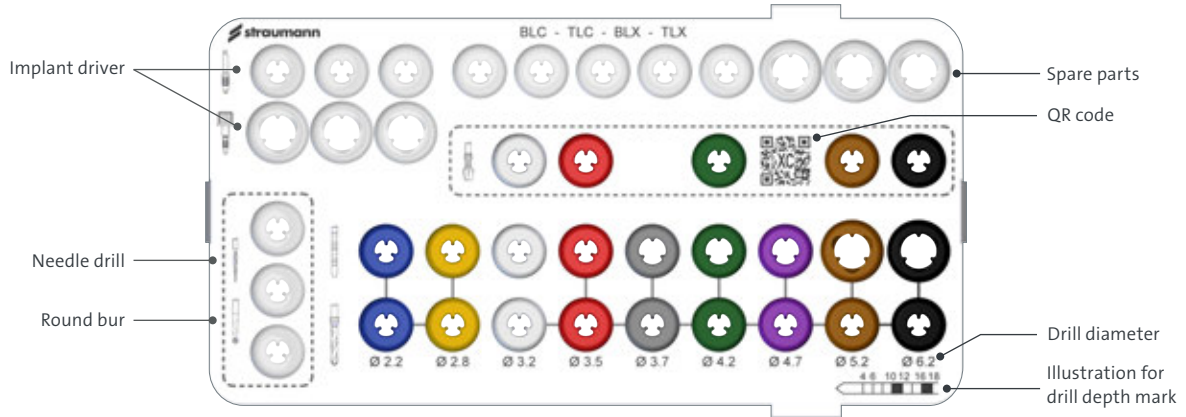
The Straumann iEXCEL™ Tray is intended for BLC, TLC, BLX, TLX implants using freehand surgery. For complete information on the surgical and prosthetic procedures for the Straumann BLC™/TLC™/BLX™/TLX™ Implant System, please refer to *Straumann iEXCEL™ Implant System, Basic Information (707743/en)*.

041.787 BLC, TLC, BLX, TLX Tray




The Straumann iEXCEL™ Tray 2 lines of drills Short & Long is intended for BLC, TLC, BLX, TLX implants using freehand surgery. For complete information on the surgical and prosthetic procedures for the Straumann BLC™/TLC™/BLX™/TLX™ Implant System, please refer to *Straumann iEXCEL™ Implant System, Basic Information (707743/en)*.

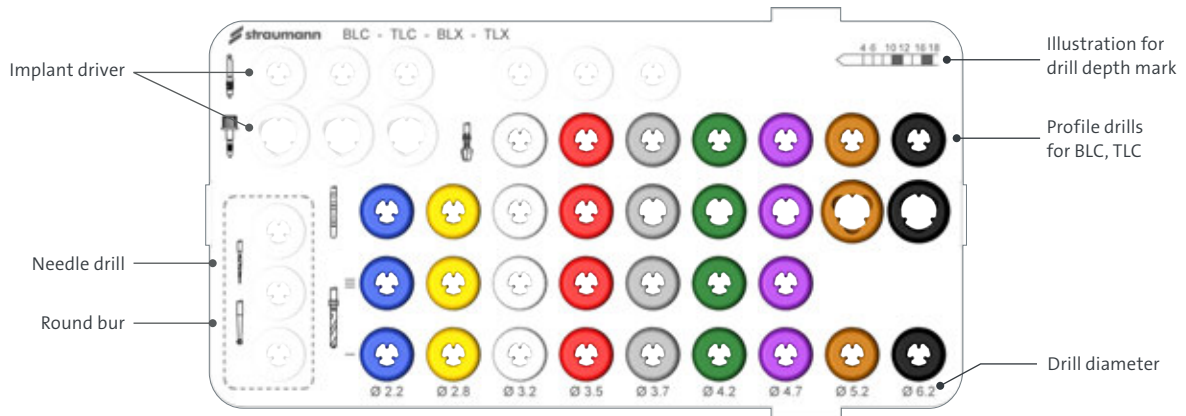041.814 BLC, TLC, BLX, TLX 2 lines of drills Short & Long Tray




The Guided Straumann iEXCEL™ Tray is intended for BLC, TLC, BLX, TLX implants using guided surgery. For complete information on the guided surgery procedure, please refer to *Straumann iEXCEL™ Guided Surgery, Basic Information (707745/en)*.

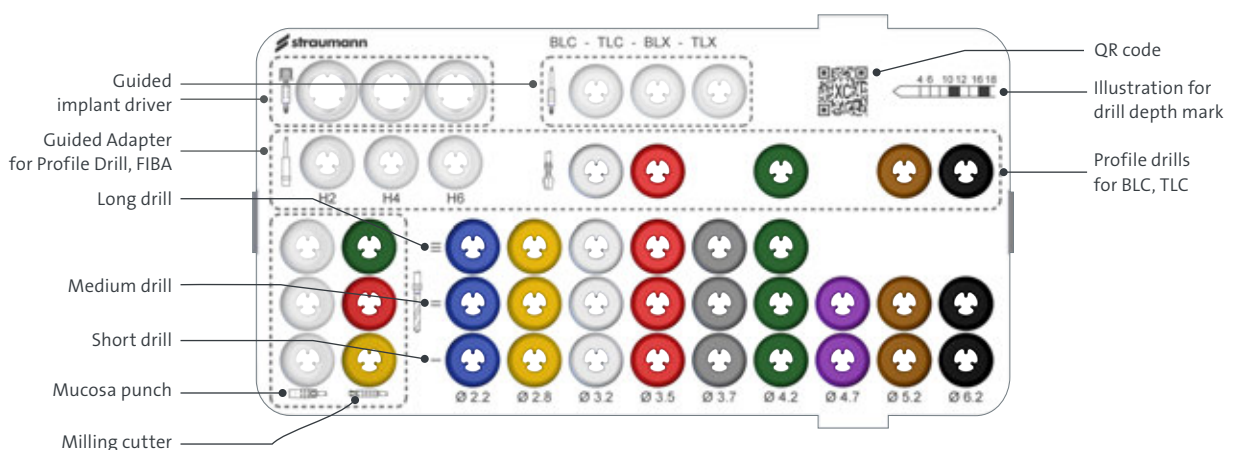041.788 Guided BLC, TLC, BLX, TLX Tray




The Straumann iGuide™ Tray is intended for BLC, TLC, BLX, TLX implants using Straumann iGuide™. For complete information on the guided surgery procedure, please refer to *Straumann iGuide™, Basic Information (707930/en)*.

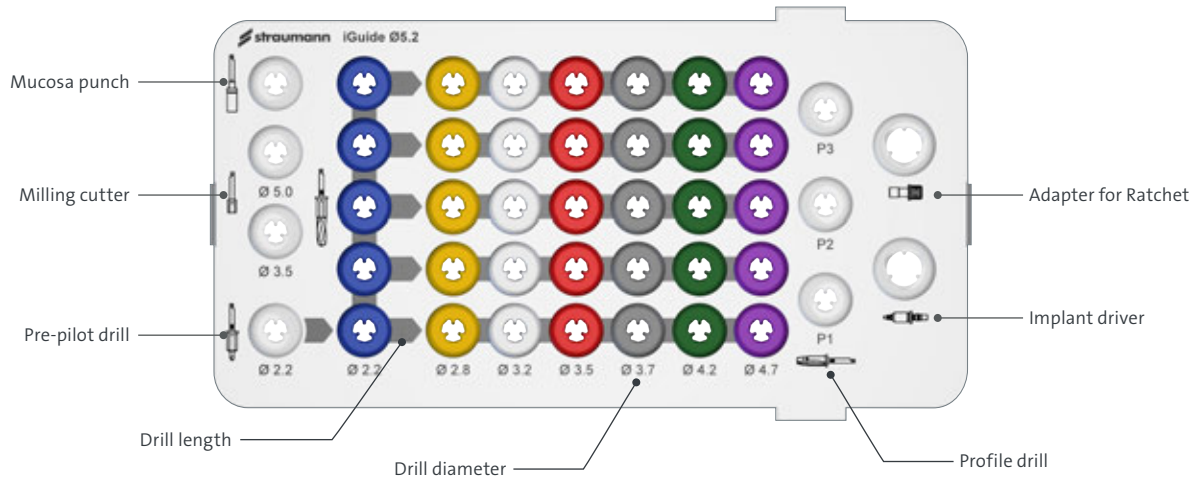041.789 Straumann iGuide™ B-Tray, Ø5.2




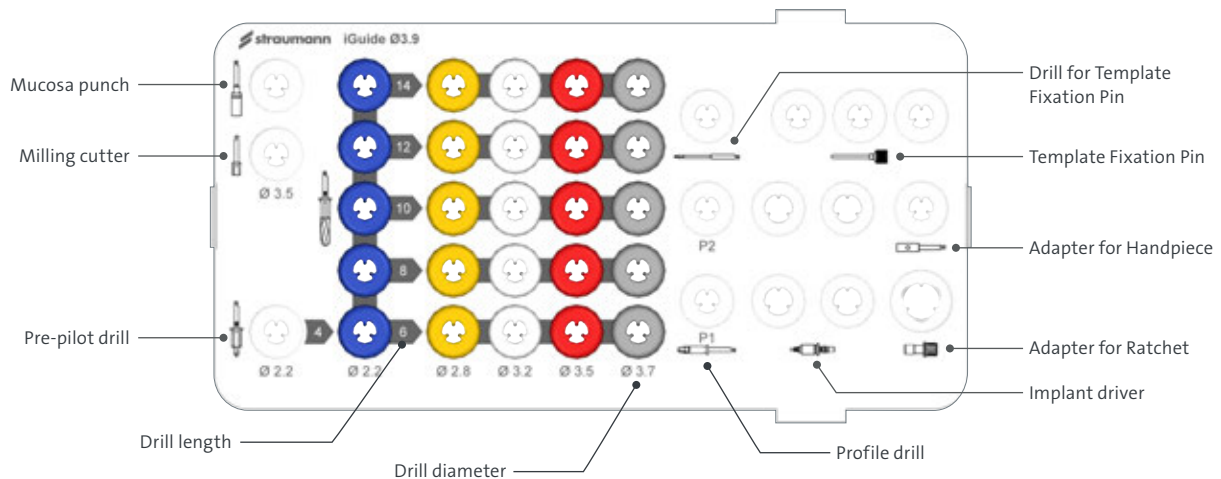041.811 Straumann iGuide™ B-Tray, Ø3.9



The BLC/TLC Tray is intended for BLC and TLC implants using freehand surgery. For complete information on the surgical and prosthetic procedures for the Straumann BLC™/TLC™ Implant System, please refer to *Straumann BLC™ and TLC™ Implant System, Basic Information ((705650 & 705651)/en)*.

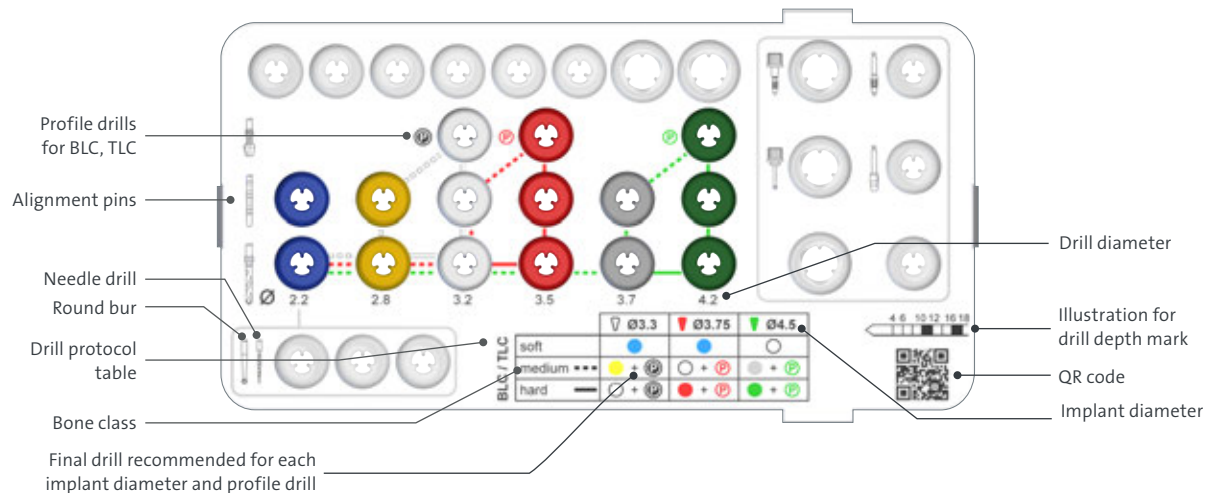041.785 BLC, TLC Tray




### 1.2.2 Fully Tapered Tray and Guided Fully Tapered Tray


The Fully Tapered Tray is intended for BLX and TLX implants using freehand surgery. For complete information on the surgical and prosthetic procedures for the Straumann BLX™ Implant System, please refer to *Straumann BLX™ Implant System, Basic Information (702115/en)* and *Straumann TLX™ Implant system, Basic information (702854/en)*.

The Guided Fully Tapered Tray is intended for BLX and TLX implants using guided surgery. For complete information on the guided surgery procedure, please refer to *Straumann® VeloDrill™ Guided Surgery, Basic Information (702526/en)*.

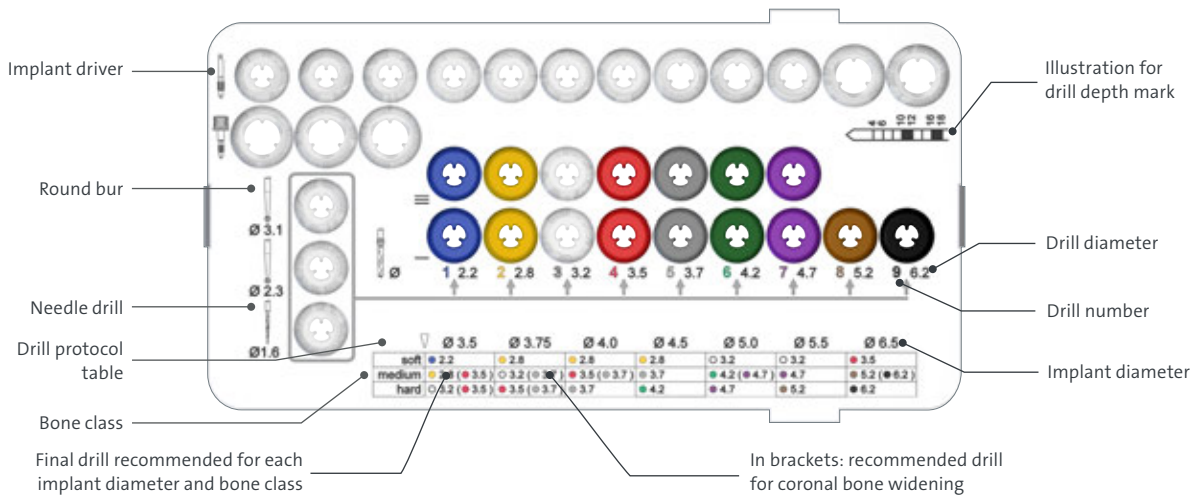041.777 B Module, Fully Tapered Tray




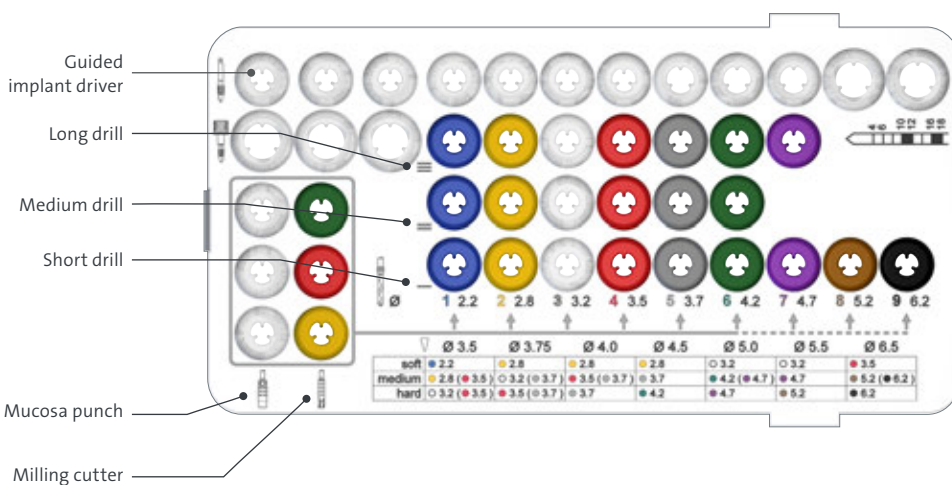
041.780 B Module, Guided Fully Tapered Tray



### 1.2.3 Basic Tray and Guided Basic Tray

The Basic Tray is intended for one implant line (e.g. SP, BL, BLT) using freehand surgery. For complete information on the surgical and prosthetic procedures, refer to *Basic information on the surgical and prosthetic procedures - Straumann® PURE Ceramic Implant (702406/en)*, *Straumann® Dental Implant System, Basic information (702084/en)* and *Straumann® Bone Level Tapered Implant, Basic Information (702167/en)*.

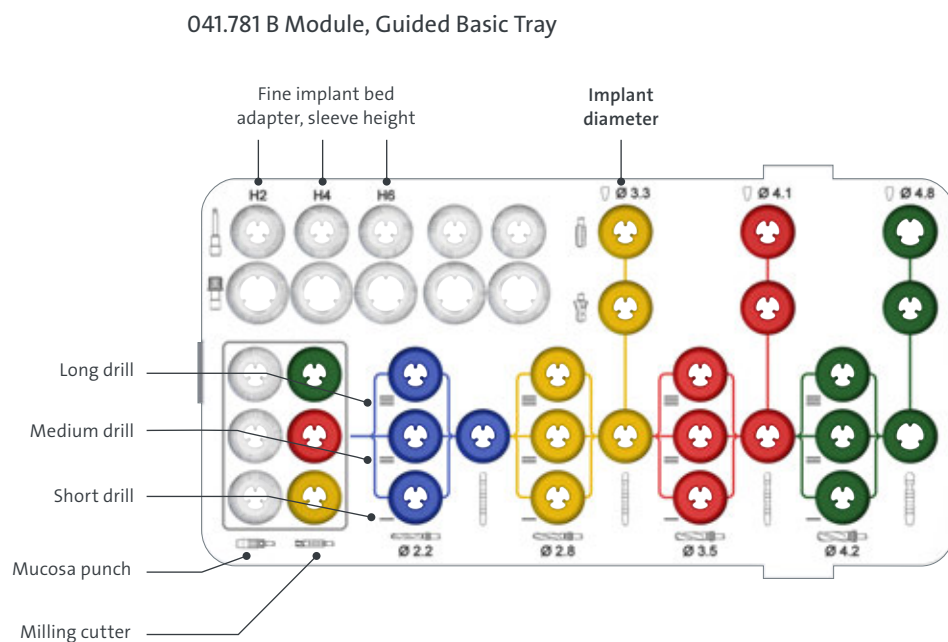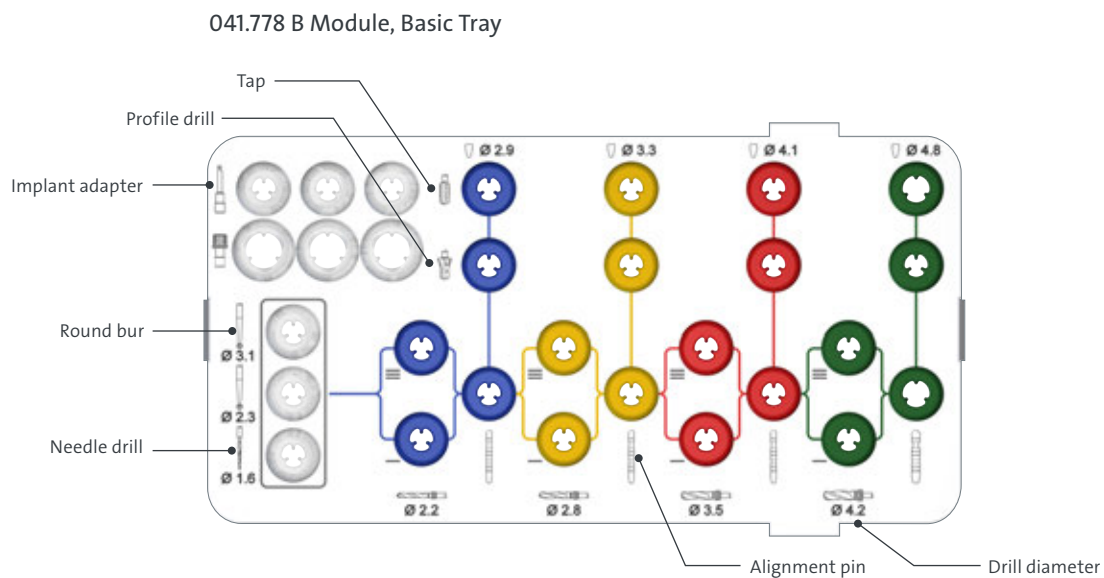The Guided Basic Tray is intended for one implant line using guided surgery. For complete information on the guided surgery procedure, please refer to *Straumann® VeloDrill™ Guided Surgery, Basic Information (702526/en)*.




### 1.2.4 Parallel Walled Tray and Guided Parallel Walled Tray


The Parallel Walled Tray is intended for Standard Plus, Bone Level and Standard Plus Short implants using freehand surgery. Refer to *Straumann® Dental Implant System, Basic information (152.754/en)* for surgical information.

The Guided Parallel Walled Tray is intended for Standard Plus, Bone Level and Standard Plus Short implants using guided surgery. For complete information on the guided surgery procedure, please refer to *Straumann® VeloDrill™ Guided Surgery, Basic Information (702526/en)*.

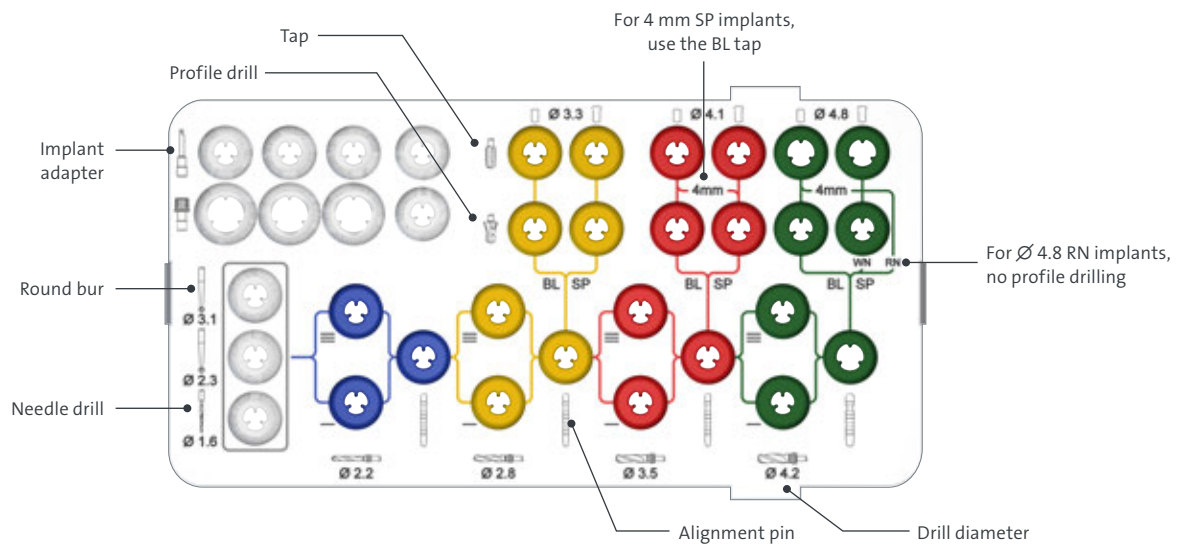041.779 B Module, Parallel Walled Tray




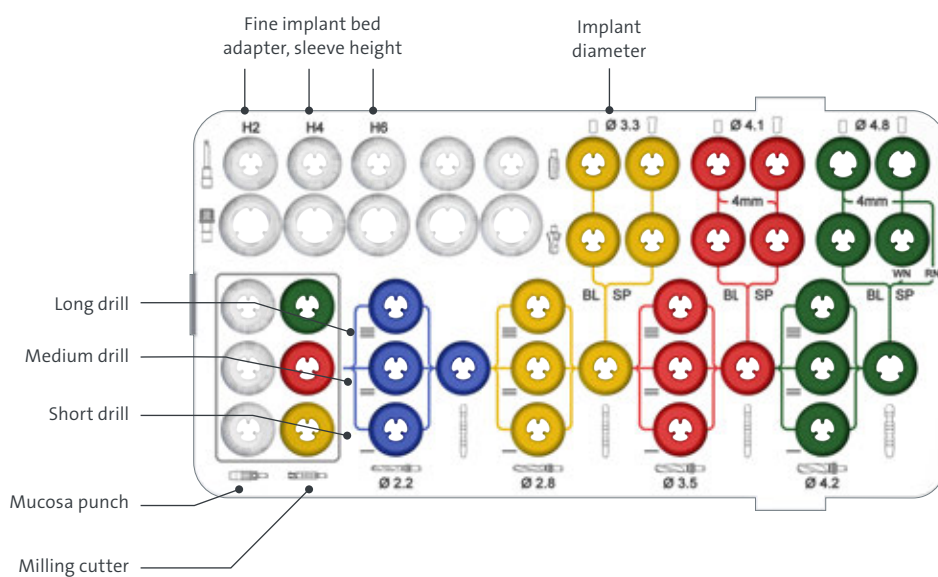041.782 B Module, Guided Parallel Walled Tray




### 1.3 STRAUMANN® MODULAR CASSETTE, C MODULE

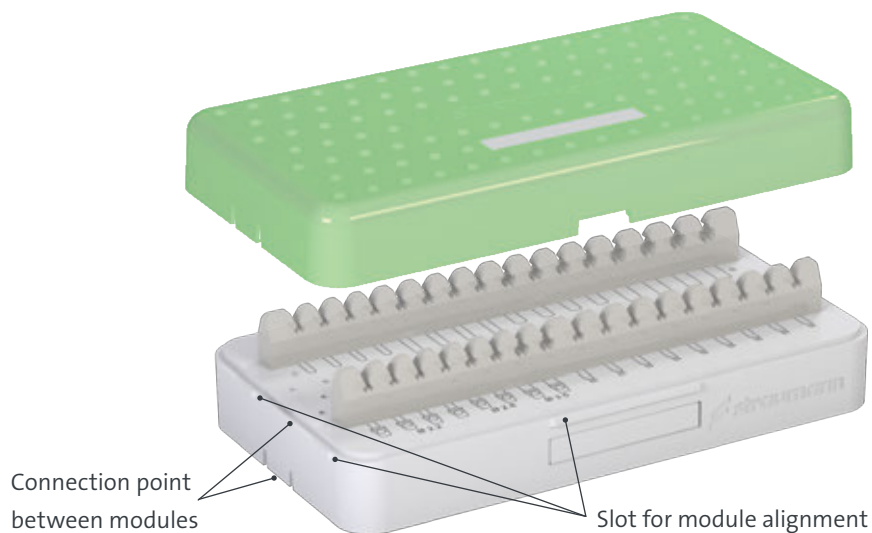C Modules have silicone strips that allow horizontal storage of instruments. External slots on the side of the module help align each module on the surgical table. Refer to the *Straumann® Modular Cassette Selection Guide (702824/en)* to select the correct cassette set-up for your clinic.



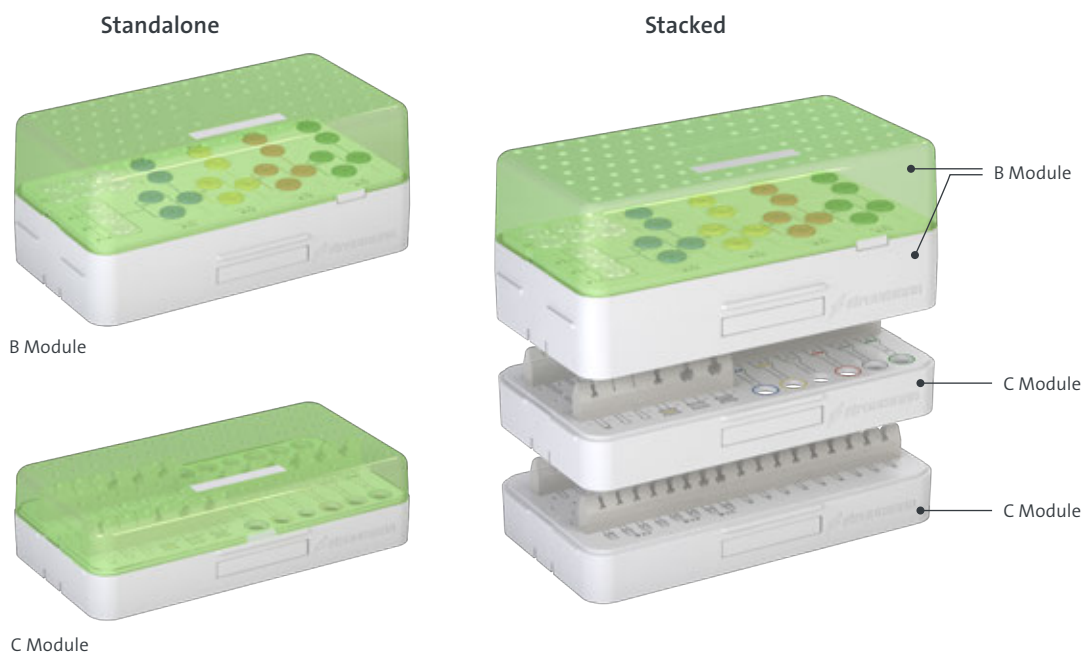The C Module can be used on its own, stacked with other C Modules or connected under a B Module base for sterilization and storage. The maximum permissible stacking height is one B Module on top of two C Module bases, without the C Module lids.



Gently push the different bases by hand to engage. Physical feedback and an audible click will confirm secure placement. To remove the base, gently pull on one side by hand to break the connection. Note that the modules are easiest to separate when the force is evenly distributed over the long side of the modules, as shown in the picture.

**Assemble**

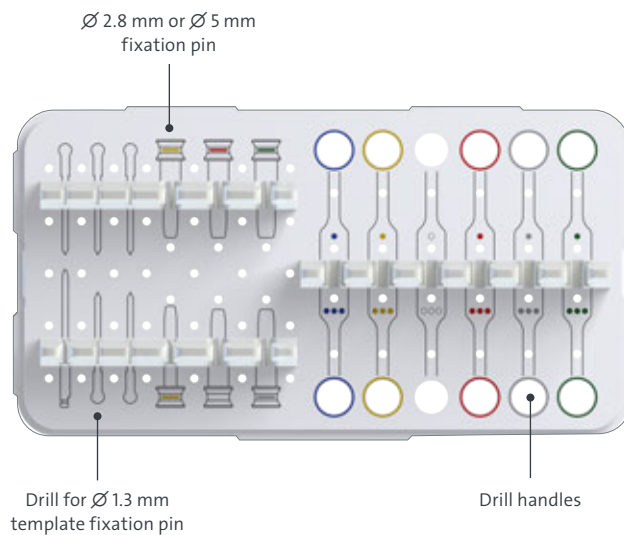


**Disassemble**



**1.3.1 C Module Guided Surgery (041.772)**


The guided surgery module stores guided handles and instrumentation for template fixation pins. All tools are stored horizontally on the holder.



### 1.3.2 C Module Spare Tools (041.773)

The spare tools module can be used to store other auxiliaries for implant surgery. For example:

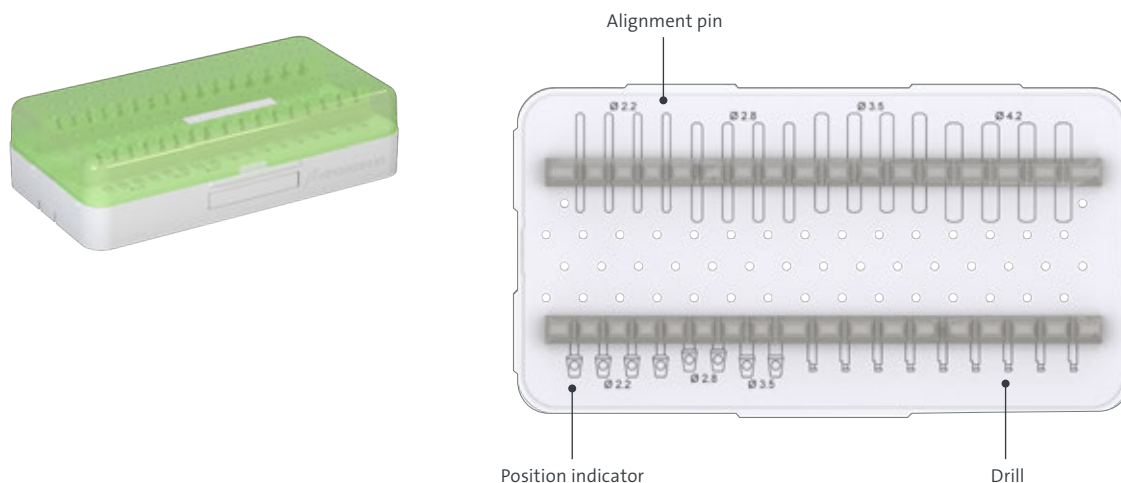- Straumann PURE™ position indicator for the Straumann® PURE Ceramic Implant
- Profile drills and taps for a specific implant line
- Additional alignment pins for Straumann ProArch™ cases.




## 1.4 LIDS

The B Module lid can be used in combination with an Ultrasonic Mat (Art. No. 041.774) for ultrasonic bath cleaning of contaminated instruments. After surgery, contaminated instruments should be pretreated and then placed on the Ultrasonic Mat that is placed in the B Module lid for ultrasonic bath cleaning.

The C Module lid can be used to enclose the B Module lid to protect the contaminated instruments from dropping out during transportation between uses and reprocessing.

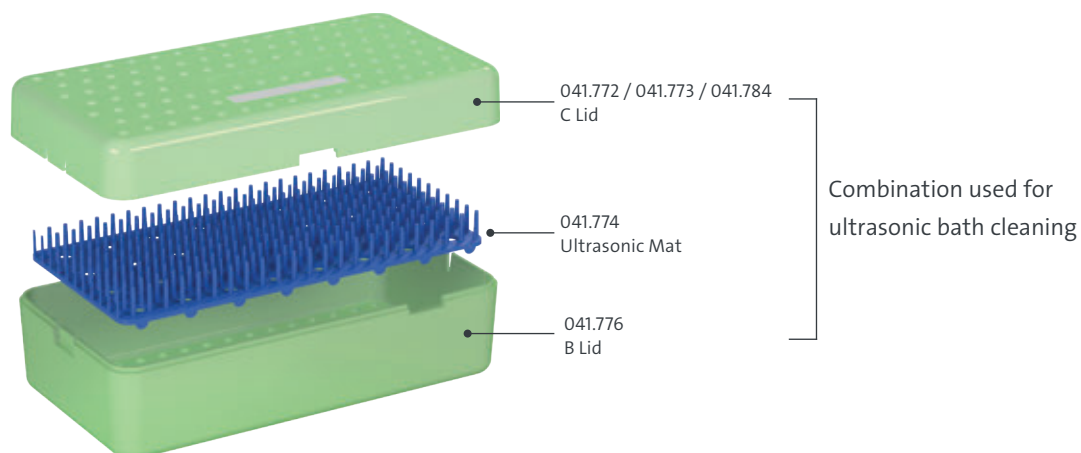For more information on the handling and reprocessing of Straumann® Instruments, please refer to chapter 2 – *Step-by-step reprocessing*, and *Straumann® Surgical and Prosthetic instruments, Care and maintenance* (152.008/en).



# 2. STEP-BY-STEP REPROCESSING

Always wear protective clothes for your own safety. Use only properly maintained equipment and materials approved according to national laws and regulations.

## 2.1 POINT OF USE



To avoid injury during surgery, the cutting instruments should be handled with the instrument tweezers (Art. No. 046.820).

Make sure that all contaminated instruments are collected separately. The instruments can be damaged by incorrect handling, such as throwing them on a hard surface. Do not reinsert the contaminated instruments back into the trays and grommets.



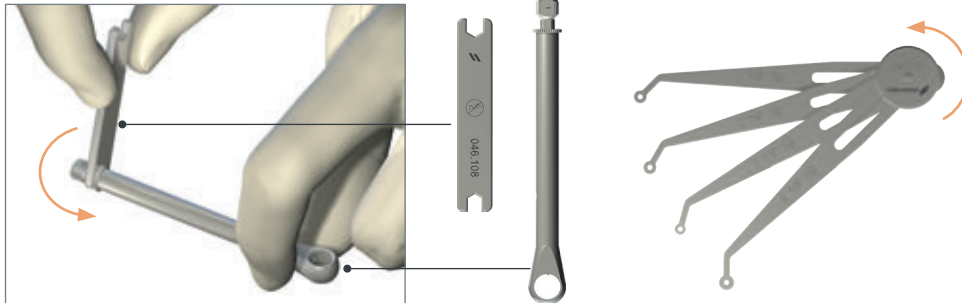
Damaged and/or blunt instruments must be sorted out and disinfected, cleaned and disposed of separately. Process contaminated instruments as soon as possible for cleaning (within one (1) hour at the most). Immerse instruments in disinfection solution after usage to avoid drying of debris.

## 2.2 CONTAINMENT AND TRANSPORTATION



Safely store and transport the cassette with instruments in a closed container to the reprocessing area to avoid any damage to, or contamination of, the environment.

## 2.3 PRETREATMENT



Sort the instruments according to their material groups to clean these groups separately.

Disassemble all trays, modules and multi-piece instruments into their individual parts (e.g. ratchet (Art. No. 046.119), distance indicator, 48h Explanation Device). To disassemble the ratchet (Art. No. 046.119), loosen the cover screw with the service instrument for ratchet (Art. No. 046.108) and unscrew the ratchet bolt.



A Module

B Module

C Module

The modular cassette should be disassembled into each individual component (lid / base / tray). To remove trays from the A Module, use a blunt instrument (e.g. the ratchet bolt from the disassembled ratchet, Art No. 046.119) to push the trays out from the back. For the B Module trays, gently pull the tray out by hand. See chapter 1 *Straumann® Modular Cassette – System overview* for more information on assembling and disassembling.



Immerse the disassembled trays, modules and multi-piece instruments in a water bath for 10 minutes.

All visible dirt should be removed by brushing with a suitable soft bristled brush under running tap water. Never use metal brushes or steel wool.

All movable components should be actuated 3 times under running water.

Use the cleaning brush for ratchet (045.111V4) or an appropriate soft nylon bottle brush to clean the lumina of the ratchet and the torque control device.

## 24 CLEANING

When selecting the cleaning agents, ensure the following:

- Always follow the instructions for use of the manufacturers of cleaning agents.
- Cleaning agents are suitable for cleaning metal and non-metallic instruments.
- The cleaning agent (if used) is suitable for ultrasonic cleaning (no development of foam).
- The chemicals used are compatible with the instruments and cassette components.
- Do not combine different cleaning agents.
- The concentrations and action times stated by the manufacturer of the cleaning agent must be strictly adhered to. Use only freshly made solutions.

Pollutants/minerals from water can reduce the lifetime of instruments and impair the performance of cleaning solutions. Always use the highest permissible water quality for cleaning (e.g. distilled or deionized). At least potable water quality should be used. For the final rinse Straumann® strongly recommends the use of fully demineralized, endotoxin-free water.

Use only filtered air for drying.

## 2.4.1 Manual cleaning and disinfection procedure



### Manual Cleaning

- Place the disassembled cassette parts (lid, trays, base) in an ultrasonic bath.
- Place the disassembled instruments on an Ultrasonic Mat (041.774) that is secured in a B Module Lid. The C Module lid can be used to enclose the B Module lid to protect the contaminated instruments during transportation.
- Make sure that the instruments do not touch each other. Instruments made of different materials may not be placed in the same bath.
- Run an ultrasonic cleaning cycle (frequency 35 kHz) for 10 minutes in a bath of deionized water supplemented with a suitable cleaning agent (e.g. 0.8% Cidezyme®, ASP).
- Rinse out all cavities of the instruments and cassette parts three times (3×) under running deionized water for 10 seconds or until no dirt residues are visible.



### Manual Disinfection

- Place the disassembled cassette parts in a disinfection bath for the specified action time (e.g. 12 min in CIDEX® OPA) at room temperature (20 °C ± 2 °C). Ensure that the cassette parts are covered by the disinfection solution.
- Rinse out all cassette parts three times (3×) with disinfection solution at the beginning and at the end of the action time using a disposable syringe (minimum syringe volume 20 ml).
- Remove the instruments from the disinfection bath.

### Rinsing

- Immerse cassette parts completely in a large volume (e.g. 8 litres) of water for a minimum of 1 minute.
- Rinse cassette parts thoroughly with water at least five times (5×). Rinse out all cavities with water and lumens of the instruments five times (5×) using a disposable syringe (minimum syringe volume 20 ml).
- Remove the cassette parts and discard the rinse water. Always use fresh volumes of water for each rinse. Do not reuse the water for rinsing or any other purpose.
- Repeat rinsing steps above for a total of THREE (3) RINSES, with large volumes of fresh water to remove CIDEX® OPA solution residues. Residues may cause serious side effects.
- Dry the instruments and cassette parts inside and outside with filtered compressed air.
- Inspect and pack the instruments and cassette parts as quickly as possible after removal. If additional drying is necessary, dry in a clean location.

## 2.4.2 Automated cleaning and disinfection (EN ISO 15883-1)



When using a washer-disinfector for automated cleaning and thermal disinfection, ensure the following:

- The washer-disinfector has been tested for its effectiveness.
- A tested program for thermal disinfection (A0 value > 3000 or) is used.
- The program used for the instruments is suitable and contains enough rinsing cycles.
- The air used for drying is filtered.
- The washer-disinfector is regularly maintained and checked.

### Procedure

- After pretreatment, mount the A module trays into the A module while leaving the A module open. Mount the B module tray into the B module without enclosing the lid.
- Place the modules in the washer-disinfector flat and opened. Do not mount lids onto the modules. Do not stack modules. Do not load instruments into the cassette. Place the ultrasonic mat alone.
- Start the appropriate program suitable for instruments. Use appropriate cleaning agents according to the manufacturer's instructions (e.g. 0.5% (v/v) neodisher® MediZym, Dr. Weigert).
- Perform thermal disinfection at 90°C for 5min.
- Remove the cassette from the washer-disinfector after the end of the program.
- Inspect and pack the cassette as quickly as possible after removal. If additional drying is necessary, dry in a clean location.

### Note

Always follow the instructions stated by the manufacturer of the washer-disinfector.

The following devices, materials and machines have been used in the cleaning and disinfection validation study:

### Manual

**Cleaning agent:** 0.8% Cidezyme®, ASP

**Disinfectant agent:** CIDEX® OPA, ASP

### Automated

**Cleaning agent:** 0.5% neodisher® MediZym, Dr. Weigert

**Washer/Disinfector:** Miele G 7836 CD

**Instrument Rack:** Miele E327 two level rack

## 2.5 INSPECTION

Check all instruments after cleaning for corrosion, damaged surfaces, chipping and contamination and sort out damaged instruments. Critical areas such as handle structures, joints or blind holes, in particular, must be inspected carefully. You can use a magnifying glass and direct lighting for better visibility. Instruments with illegible markings/labeling must also be replaced.

Instruments which are still contaminated must be cleaned again. Damaged, corroded or worn instruments should not come into contact with intact instruments to avoid contact corrosion.

Check the modular cassette for functionality:

- Inspect the condition of all modules. Replace modules if there is damage, markings, deformation or wear.
- Visually inspect that all labeling printed on the lid, tray and modules is clear and legible
- Ensure that all trays are correctly assembled with the cassette before loading with instruments. The trays are connected correctly when they do not fall out of the cassette when the cassette is turned upside down.
- Ensure a secure connection between lid and module base.
- For stacked storage, ensure a secure connection between modules
  - B Module base on C Module base
  - C Module base on C Module base

## 2.6 PACKAGING



Make sure that the cassette and instruments are completely dry before packing for sterilization.

Place all instruments and the disassembled ratchet (Art. No. 046.119) in the designated slots in the cassette. Assemble the cassette by putting together the tray, module and lid. Refer to the *Straumann® Modular Cassette Selection Guide* (702824/en) for guidance on loading the cassette.

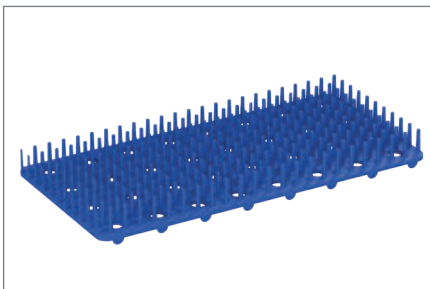


When connecting different modules, the maximum permissible stacking limit is one B Module on top of two C Module bases without C lids. The A Module should be sterilized alone.

An indicator strip with the date of the sterilization should be affixed to every sterilization packaging. This will help to indicate whether and, if so, when the material was sterilized.

Place the cassette in a metal sterilization container or in double-pouch packaging satisfying the following requirements:

- Suitable for steam sterilization
- EN ISO/ANSI AAMI ISO 11607 – Packaging for terminally sterilized medical devices
- For USA: Use FDA-cleared sterilization accessories



The Ultrasonic Mat should be packaged and sterilized individually in a double pouch. Sterilizing instruments on the Ultrasonic Mat is not allowed.

## 2.7 STERILIZATION

As the cassette is not intended to maintain sterility on its own, place the instruments in the appropriate module and seal the module in a sterilization double pouch or metal sterilization container.

When loading the sterilizer, always place the cassette on the shelf in such a way that it does not come in contact with the walls of the sterilizer. Do not put the cassette on its side or upside down on top of the lid.

Do not place corroded, rusty instruments in the cassette for sterilization. These contaminate the water circulation system of the sterilizer with rust particles. During every subsequent sterilization cycle, these rust particles cause rust on instruments that were originally intact.

Only steam sterilization methods listed below may be used for sterilization. Other sterilization methods are not allowed. The sterilizer manufacturer's instructions for use must be strictly followed.

The sterilizer must:

- Offer a fractionated vacuum method with sufficient device drying time and compliant with EN 13060 or EN 285
- Be validated according to EN ISO 17665 (valid IQ/OQ and product-specific performance assessment (PQ))
- Provide a maximum sterilization temperature of 134 °C (273 °F; plus tolerance corresponding to DIN EN ISO 17665, i.e. 137 °C (278 °F)).

### Recommended sterilization time (exposure time at the sterilization temperature) and drying time:

Method	Conditions	Drying time
<b>For Europe:</b>		
Moist Heat (Autoclave) Fractionated vacuum	134 °C (273 °F) for 3 min	Local practice
<b>For the United States:</b>		
Moist Heat (Autoclave) Fractionated vacuum	132 °C (270 °F) for 4 min	30 min
<b>For countries outside Europe and the United States:</b>		
Moist Heat (Autoclave) Fractionated vacuum	132 °C (270 °F) to 134 °C (273 °F) at least for 3 min	Local practice

If visible signs of moisture are present (damp spots on sterile packaging, pooled water in the load) at the end of the sterilization cycle, repackage and re-sterilize using a longer drying time.

**Note:** Users in the United States should ensure that the sterilizer and all sterilization accessories (sterilization wraps, pouches, sterilization trays, biological indicators and chemical indicators) are cleared by the FDA for the intended sterilization cycle.

# 3. FURTHER INFORMATION






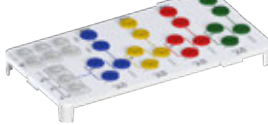

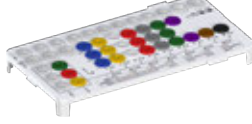





For further information (e.g. warnings, cautions, precautions, compatibility), please consult the *Instructions for use Straumann® Modular Cassette (702407/en)*, *Straumann® Modular Cassette Selection Guide (702824/en)* and *Straumann® Surgical and Prosthetic Instruments, Care and maintenance (702000/en)*.

# 4. PRODUCT REFERENCE LIST



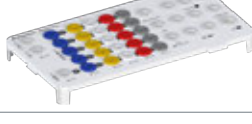
## A Module overview

A Module overview		
Art. No.	Picture	Description
041.761		Straumann® Modular Cassette, A Module
041.762		A Module, Grommet Tray, 6 small
041.763		A Module, Grommet Tray, 6 large
041.764		A Module, Grommet Tray, 3 small 3 large
041.765		A Module, Grommet Tray, 3 medium 3 large
041.791		A Module, Ratchet with Torque Control Tray
041.766		A Module, Ratchet Tray
041.767		A Module, Release Aid Tray
041.768		A Module, Stop Key Tray
041.769		A Module, Planning Aid Tray
041.770		A Module, Tweezer Tray
041.771		A Module, Implant Depth Gauge Tray

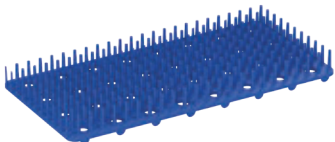

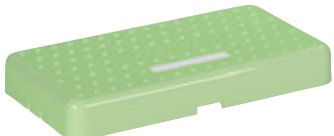
## B Module overview

B Module overview		
Art. No.	Picture	Description
041.776		Straumann® Modular Cassette, B Module
041.785		B Module, BLC, TLC Tray
041.777		B Module, Fully Tapered Tray
041.814		B Module, BLC, TLC, BLX, TLX, Short & Long Drills Tray
041.787		B Module BLC, TLC, BLX, TLX Tray
041.778		B Module, Basic Tray
041.779		B Module, Parallel Walled Tray
041.780		B Module, Guided Fully Tapered Tray
041.788		B Module, Guided BLC, TLC, BLX, TLX Tray
041.781		B Module, Guided Basic Tray
041.782		B Module, Guided Parallel Walled Tray
041.789		Straumann iGuide™ B-Tray, Ø5.2
041.811		Straumann iGuide™ B-Tray, Ø3.9

## C Module overview

C Module overview		
Art. No.	Picture	Description
041.783		Straumann® Modular Cassette, C Module, Auxiliary Box
041.772		Straumann® Modular Cassette, C Module Guided surgery
041.773		Straumann® Modular Cassette, C Module Spare Tools

## Auxiliaries

Auxiliaries		
Art. No.	Picture	Description
041.774		Ultrasonic Mat
041.775		Spare Label Plate: Blank BLX BLT SP BL BLC/BLX/TLC/TLX SP/BL Short Implant BLC/TLC Auxiliary Box BLX Guided BLT Guided SP Guided BL Guided BLC/BLX/TLC/TLX Guided SP/BL Guided Short Implant Guided Short Implant Guided Guided Surgery
041.784		C Module Lid







## **International Headquarters**

Institut Straumann AG

Peter Merian-Weg 12

CH-4002 Basel, Switzerland

Phone +41 (0)61 965 11 11

Fax +41 (0)61 965 11 01

[www.straumann.com](http://www.straumann.com)

© Institut Straumann AG, 2026. All rights reserved.

Straumann® and/or other trademarks and logos from Straumann® mentioned herein are the trademarks or registered trademarks of Straumann Holding AG and/or its affiliates.

